

DONJOY®

IsoFORM™

CUSTOMIZABLE SPINE
BRACES DESIGNED
TO BE AS UNIQUE AS
YOUR PATIENTS.



enovis™




DONJOY[®]

IsoFORM™ Back Braces

Comfortable compression for lower back stability and support

Designed to aid superior lower back protection and stabilization, DonJoy[®] IsoFORM™ back braces help provide support and pain relief for a comprehensive range of chronic and acute indications. Supplement support with hot/cold packs for additional relief.

Built with customizable support, targeted compression, and patient comfort in mind

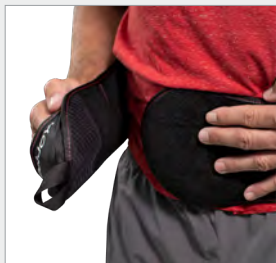


Adjustable Belt Wings

- accommodate different shapes for a custom fit



Paneling

- provides durable rigid support

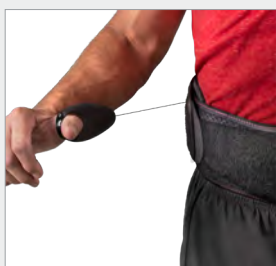


Dual Pulley System

- provides targeted compression with independent top and bottom compression



Flexible Pull Tabs

- for ease of use with adjustable lace length feature

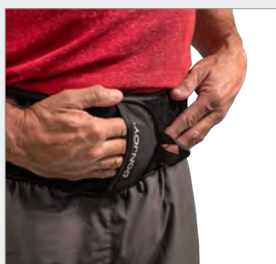

Donning Mitt

- assists with ease of application



Semi-Universal

- comfortably fits waist sizes from 71 - 152 cm




Hexagonal Spacer Mesh

- durable and breathable materials for comfort




Optional Hot/Cold Pack

- for additional relief

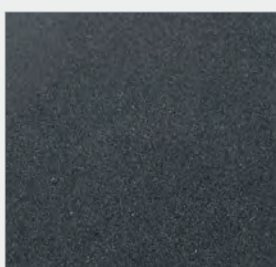

EVA-lined Panels

- provides comfortable support

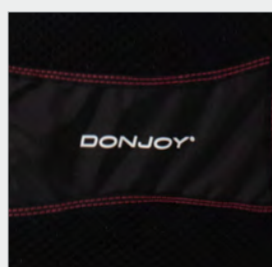

Enclosed Lacing panel

- protects lacing system from unwanted snags or tears



IsoFORM® SIO

- Flexible, low profile belt with single pulley system provides compression to reduce SI joint pain and dysfunction. Designed with silicone grip banding to help prevent brace from migrating.

Indications: Suitable for people with sacroiliac pain, dysfunction, lower back pain, sprains/strains

IsoFORM – LO/LO+

IsoFORM LO

- Provides sagittal control and support from L1-L5. With optional Hot Cold Pack for additional relief.

IsoFORM LO+

- Similar design as the LO, plus a Removable Semi-Rigid Anterior Panel for additional trunk support.

Indications: Designed to support and stabilize the lower back and may be suitable for people with; acute and chronic low back pain, lower back sprains or strains, osteoporosis, spondylolisthesis



IsoFORM – LSO/LSO+

IsoFORM LSO

- With posterior panel and removable anterior panel provides sagittal control and support from T9-S1.

IsoFORM LSO+

- Similar design to the LSO, plus two removable semi-rigid lateral extensions to help reduce lateral movement.

Indications: Designed to support and stabilize the lower back and may be suitable for people with acute and chronic low back pain, lower back sprains or strains, osteoporosis, spondylolisthesis



IsoFORM® Postural Extension



Customizable Support

- Telescoping aluminum posterior support and shoulder straps are adjustable and malleable for custom fit
- Removable semi-rigid anterior panel for additional trunk support
- Sternal cross strap secures the shoulder straps for a comfortable fit
- Adjustable belt wings accommodate different shapes for a custom fit

Targeted Compression

- Dual pulley system provides independent top and bottom compression
- Enclosed lacing panel protects lacing system from unwanted snags or tears

Patient Comfort

- Adjustable backpack-style straps for easy donning and doffing
- Donning mitt pocket assists with application
- Flexible pull tabs – for ease of use with adjustable lace-length feature
- Optional hot/cold pack for additional relief

Indications: Suitable for people with thoracic spine conditions such as kyphosis, osteoporosis, acute pain, mechanical back pain, compression fractures, strain and mild scoliosis



A. Adjustable Straps Backpack style for easy donning and doffing

B. Removable Semi-Rigid Anterior Panel For additional trunk support

C. Sternal Cross-Strap For a comfortable fit

D. Adjustable Belt Wings Accommodate different shapes and allow clinicians to customize the fit for a wide range of patients



E. Telescoping Aluminum Posterior Support Height adjustable and malleable to help tailor to the patient

F. Enclosed Lacing Panel Protects lacing system from unwanted snags or tears

G. Dual Pulley System With independent top and bottom compression

H. Pull Tabs For comfortable donning and doffing

IsoFORM® TLSO

Customizable Support

- Telescoping aluminum sternal bar, adjustable and malleable for custom fit
- Anterior panel helps restrict flexion and trunk rotation and supports T3-S1
- Posterior panel – rigid plastic helps restrict extension and trunk rotation and supports T3-S1
- Lateral extensions - removable semi-rigid supports help reduce lateral movement
- Belt wings with removable plastic paneling encompass the torso for customizable support
- Accommodate different shapes for a custom fit

Targeted Compression

- Dual pulley system provides independent top and bottom compression
- Enclosed lacing panel protects lacing system from unwanted snags or tears

Patient Comfort

- Easy to apply and remove with pivoting quick-release buckles
- Over-the-shoulder or underarm strap configuration options
- Donning mitt pocket assists with application
- Flexible pull tabs for ease of use with adjustable lace-length feature
- Optional hot/cold pack for additional relief




Indications: suitable for people with compression fractures of the thoracic spine, osteoporosis, thoracic mechanical back pain and kyphosis



- A. Pivoting Quick Release Buckles** For easy donning and doffing
- B. Aluminum Sternal Bar** Telescoping bar can be adjusted to desired height and malleable for custom fitting to desired angle
- C. Semi-Rigid Belt Wings** Encompass the torso for added support (removable plastic panel for customizable support)
- D. Lateral Extensions** Semi-rigid supports to help reduce lateral movement (Removable for customizable support)
- E. Anterior Panel** Helps restrict forward flexion and trunk rotation



- F. 2 Strapping Configurations** Optional over-the-shoulder or underarm strap configuration
- G. Dual Pulley System** With independent top and bottom compression
- H. Enclosed Lacing Panel** Protects lacing system from unwanted snags or tears
- I. Pull Tabs** For comfortable donning and doffing
- J. Posterior Panel** Plastic posterior panel helps restrict motion

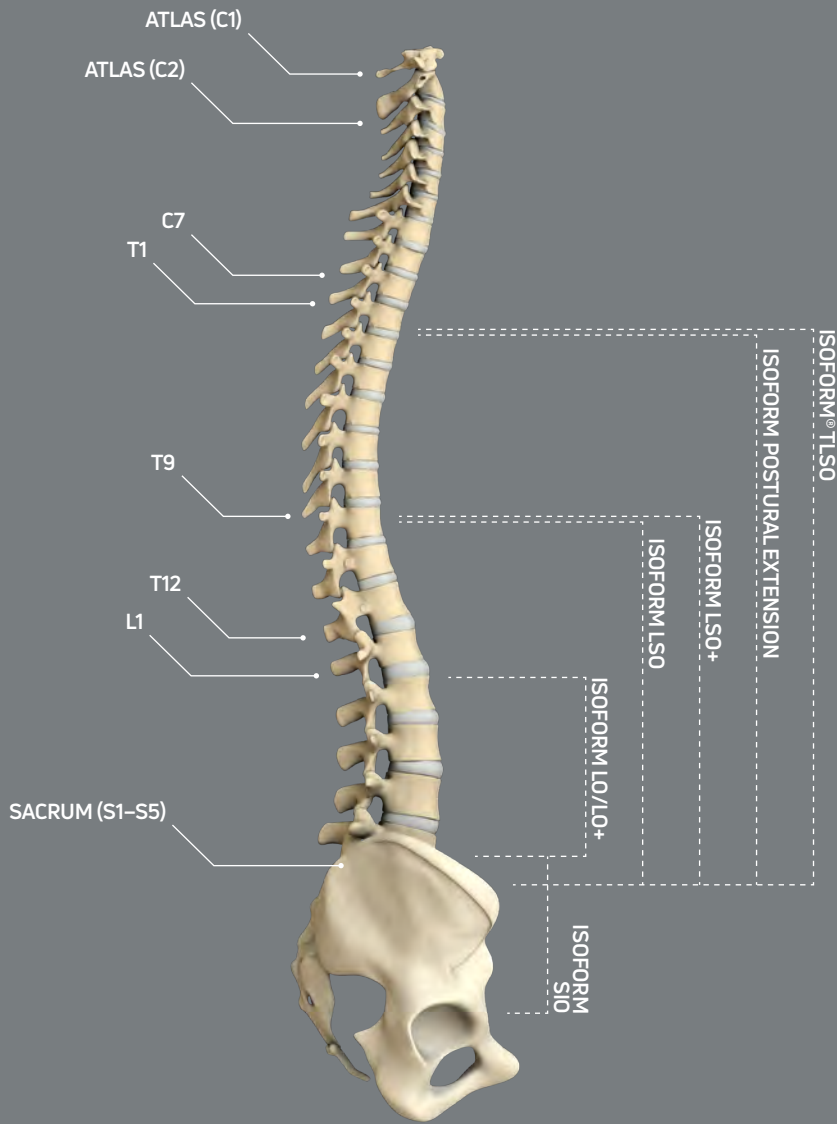

Anatomical Range of Spinal Support

DESCRIPTION	PART NUMBER	SIZE	HIP CIRCUMFERENCE
IsoFORM SIO	11-1692-2	S/M	71-122 cm (28-48")
	11-1692-4	L/XL	122-152 cm (48-60")
IsoFORM - LO	11-1693-2	S/M	71-122 cm (28-48")
	11-1693-4	L/XL	122-152 cm (48-60")
IsoFORM - LO+	11-1694-2	S/M	71-122 cm (28-48")
	11-1694-4	L/XL	122-152 cm (48-60")
IsoFORM - LSO	11-1695-2	S/M	71-122 cm (28-48")
	11-1695-4	L/XL	122-152 cm (48-60")
IsoFORM - LSO+	11-1696-2	S/M	71-122 cm (28-48")
	11-1696-4	L/XL	122-152 cm (48-60")
IsoFORM Postural Extension	11-1697-2	S/M	71-122 cm (28-48")
	11-1697-4	L/XL	122-152 cm (48-60")
IsoFORM TLSO	11-1698-2	S/M	71-122 cm (28-48")
	11-1698-4	L/XL	122-152 cm (48-60")
Sternal Support Accessory	11-1700	XS/S	N/A
	11-1701	XL/XXL	N/A
Hot/Cold Pack Accessory	11-1699	Can be used on LO through TLSO product offerings	N/A



Asia Pacific

1905, Tower 2, Grand Central Plaza
Shatin
HONG KONG
T: 852 31051415
E: marketing.ap@djoglobal.com

Australia – New Zealand

ANZ Asia
PO Box 6057
Frenchs Forest DC
NSW 2086
AUSTRALIA
T: 1300 66 77 30
E: customerservice.au@djoglobal.com

Benelux

Kleinhoefstraat 5, bus 39
Geel 2440
BELGIUM
T: +32 (0)14 248090
E: benelux.orders@djoglobal.com

Canada

6485 Kennedy Rd
Mississauga
Ontario, L5T2W4
CANADA
T: 1-866-866-5031
E: canada.orders@djoglobal.com

China

A312 SOHO ZhongShan Plaza
1055 W. ZhongShan Rd
Shanghai, 200051
CHINA
T: (8621) 6031 9989
E: information_china@djoglobal.com

France

3 rue de Bethar
Centre Européen de Frêt
64990 Mouguerre
FRANCE
T: +33 (0)5 59 52 86 90
E: sce.cial@djoglobal.com

Germany

Bötzingen Straße 90
79111 Freiburg
GERMANY
T: +49 (0) 180 1 676 333
E: kundenservice@djoglobal.de

Italy

Via Leonardo da Vinci 97
Trezzano sul Naviglio
Milan 20090
ITALY
T: +39 02 835 98001
E: vendite@djoglobal.com

Nordic

Murmansgatan 126
212 25 Malmö
SWEDEN
T: +46 40 39 40 00
E: info.nordic@djoglobal.com

South Africa

Unit 6, 2 on London, 2 London Circle
BrackenGate Business Park
Brackenfell, 7560
SOUTH AFRICA
T: +27 (0) 21 276 2968
E: info.southafrica@djoglobal.com

Spain

Av.Cornella 144. 1-4
08950 Esplugues de Llobregat
Barcelona
SPAIN
T: +34 93 480 32 02
E: ventas@djoglobal.com

UK & Ireland

1a Guildford Business Park
Guildford, Surrey
GU2 8XG
UK
T: +44 (0)1483 459 659
E: ukorders@djoglobal.com

USA

2900 Lake Vista Drive
Lewisville
TX 75067
U.S.A.
T: 760 727 1280
E: customercare@djoglobal.com

enovis[™]

Enovis
1a Guildford Business Park, Guildford, Surrey, GU2 8XG, UNITED KINGDOM

djoglobal.eu

Copyright © 2022 DJO, LLC, a subsidiary of Enovis Corporation
MKT-000-00011131-EN - Rev C

Individual results may vary. Neither Enovis, DJO, LLC or any of their subsidiaries dispense medical advice. The contents of this document do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.