

DONJOY®

Hand & Wrist range

Introducing a comprehensive family of functional hand & wrist braces for every indication



Stay in **touch.**





Hand & Wrist range

Introducing a comprehensive family of functional hand & wrist braces for every indication

Our hands form a connection to the world around us. With them we touch, hold, and make ourselves understood. Without them we feel lost.

The hand is one of the most complex parts of the body, and one of the most active. With so many bones, ligaments, tendons, and joints working together, so often, the risk of sudden or gradual injuries is considerable. This is why it is important for practitioners and patients to have a go-to solution for the most common hand and wrist injuries they might encounter.

DonJoy® presents a comprehensive family of braces designed to help patients with wrist and hand injuries to maintain the life they love. Combining ergonomic design with breathable materials, these braces provide support and comfort for patients with a wide range of indications, allowing them to continue their daily activities while they heal.

A Design for Life

This range of braces shares a variety of functional features to aid and protect the user.

Long-term protection

The use of highly breathable, skin-friendly materials enhances patient comfort and compliance for around-the-clock use



Modular concept

Removable, adjustable components allow the braces to be adapted to the patient's needs throughout the healing process

Lightweight

Thanks to the materials used, all of the braces are extremely lightweight for improved patient comfort



Tailored support

Form-fitting, adjustable designs enable the user to tailor the braces to their own needs without compromising support



Easy to fit

DonJoy's 'Stop Strap' system allows hassle-free and secure fitting with the use of only one hand

Latex free



DonJoy DuoForm®

Universal wrist immobilizer

Indications

- Wrist sprains and trauma
- Pre- and post-operative immobilization

Characteristics

- Easy application
- Universal model, fits either right or left hand
- Repositionable dorsal stay, provides appropriate tailored support
- Ergonomic palmar splint to control wrist flexion
- Removable modular components for universal use
- Wide opening for easy application
- Three straps for adjustable compression
- Breathable 3D material
- Comfortable, flexible, lightweight design

Ordering Information

Part Number	Wrist Circumference		Size
	CM	Inches	
82-96130	11 - 14 cm	4 ¼" - 5 ½"	Pediatric/XS
82-96132	14 - 16 cm	5 ½" - 6 ½"	S
82-96133	16 - 19 cm	6 ½" - 7 ½"	M
82-96134	19 - 21 cm	7 ½" - 8 ½"	L
82-96135	21 - 24 cm	8 ½" - 9 ½"	XL



DonJoy DuoForm+®

Universal wrist and thumb immobilizer

Indications

- Wrist & thumb sprains and trauma
- Pre- and post-operative immobilization

Characteristics

- Easy application
- Universal model, fits either right or left hand
- Repositionable dorsal stay, provides appropriate tailored support
- Ergonomic palmar splint to control wrist flexion
- Removable modular components for universal use
- Wide opening for easy application
- Three straps for adjustable compression
- Breathable 3D material
- Comfortable, flexible, lightweight design
- Distal thumb securing strap
- Hook and loop system for easy and secure fitting

Ordering Information

Part Number	Wrist Circumference		Size
	CM	Inches	
82-96140	11 - 14 cm	4 ¼" - 5 ½"	Pediatric/XS
82-96142	14 - 16 cm	5 ½" - 6 ½"	S
82-96143	16 - 19 cm	6 ½" - 7 ½"	M
82-96144	19 - 21 cm	7 ½" - 8 ½"	L
82-96145	21 - 24 cm	8 ½" - 9 ½"	XL



DonJoy RespiForm®

Breathable lightweight wrist brace

Indications

- Wrist sprains and trauma
- Tendinopathy
- Pre- and post-operative immobilization

Characteristics

- Dedicated for either right or left hand use
- Ergonomic, flexible lightweight design
- Anatomically designed, provides tailored individual fit
- Easy and simple application via adjustable straps and "Stop Strap" system
- Incorporates two semi-rigid dorsal stays and an aluminum palmar stay for overall wrist immobilisation
- Advanced Ozone™ material, promotes air flow
- Breathable in general use, providing good skin care and around-the-clock comfort in different heat and humidity conditions

Ordering Information

Part Number		Wrist Circumference		Size
Right	Left	CM	Inches	
82-96160	82-96150	11 - 14 cm	4 ¼" - 5 ½"	Pediatric/XS
82-96162	82-96152	14 - 16 cm	5 ½" - 6 ½"	S
82-96163	82-96153	16 - 19 cm	6 ½" - 7 ½"	M
82-96164	82-96154	19 - 21 cm	7 ½" - 8 ½"	L
82-96165	82-96155	21 - 24 cm	8 ½" - 9 ½"	XL



DonJoy RespiForm+®

Breathable lightweight wrist and thumb brace

Indications

- Wrist & thumb sprains and trauma
- Tendinopathy, rhizarthrosis
- Pre- and post-operative immobilization

Characteristics

- Dedicated for either right or left hand use
- Ergonomic, flexible lightweight design
- Anatomically designed, provides tailored individual fit
- Easy and simple application via adjustable straps and "Stop Strap" system
- Incorporates two semi-rigid dorsal stays and an aluminum palmar stay for overall wrist immobilisation
- Advanced Ozone™ material, promotes air flow
- Breathable in general use, providing good skin care and around-the-clock comfort in different heat and humidity conditions
- Modular designed malleable aluminum thumb stay for optimum positioning
- Hook and loop system for easy and secure distal thumb fitting

Ordering Information

Part Number		Wrist Circumference		Size
Right	Left	CM	Inches	
82-96180	82-96170	11 - 14 cm	4 ¼" - 5 ½"	Pediatric/XS
82-96182	82-96172	14 - 16 cm	5 ½" - 6 ½"	S
82-96183	82-96173	16 - 19 cm	6 ½" - 7 ½"	M
82-96184	82-96174	19 - 21 cm	7 ½" - 8 ½"	L
82-96185	82-96175	21 - 24 cm	8 ½" - 9 ½"	XL



DonJoy DigiForm®

Finger splint with removable stay for immobilization of 2 fingers


Indications

- All indications requiring immobilization of two fingers
- Sprains, subluxations, dislocations, fractures
- Pre- / post-surgery wrist immobilization

Characteristics

- Easy application via 4 adjustable straps and "Stop Strap" system
- Fits either right or left hand
- Lightweight, low profile design
- Prepositional palmar stay provides tailored support
- Aluminum stay is malleable to adapt to healing protocol
- Removable modular components for ease of washing
- Four removable and adjustable straps provide adjustable compression
- Finger tips protective cover
- Available in three sizes

Ordering Information

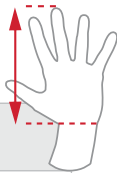

Part Number	Wrist Circumference		Size
	CM	Inches	
82-96200	11 - 15 cm	4½" - 6"	Pediatric/XS
82-96202	15 - 19 cm	6" - 7½"	S/M
82-96203	19 - 23 cm	7½" - 9"	L/XL



ANY 2
ADJACENT
FINGERS



DonJoy DigiForm+®

Finger splint with removable palmar stay for immobilization of 3-4 fingers


Indications

- All indications requiring immobilization of three or four fingers
- Sprains, subluxations, dislocations, fractures
- Pre- / post-surgery wrist immobilization

Characteristics

- Easy application via four adjustable straps and "Stop Strap" system
- Fits either right or left hand
- Lightweight, low profile design
- Prepositional aluminum palmar stay provides tailored support
- Removable modular components for ease of washing
- Finger tips protective cover
- Splint adjustable for combinations of three fingers, or for all four fingers
- Breathable 3D material
- Available in five sizes

Ordering Information



Part Number	Wrist Circumference		Size
	CM	Inches	
82-96210	11 - 14 cm	4¼" - 5½"	Pediatric/XS
82-96212	14 - 16 cm	5½" - 6½"	S
82-96213	16 - 19 cm	6½" - 7½"	M
82-96214	19 - 21 cm	7½" - 8½"	L
82-96215	21 - 24 cm	8½" - 9½"	XL



ANY 3
ADJACENT
OR ALL
FINGERS



DonJoy ErgoForm®

Universal thumb immobilizer

Indications


- All indications requiring immobilization of the thumb joints
- Indicated for sprains of the thumb, subluxations, dislocations, rhizarthrosis

Characteristics

- Ergonomic design
- Universal model, fits either right or left hand
- Easy application via adjustable strap and "Stop Strap" system
- Malleable aluminum thumb stay for optimum positioning
- Distal strap to secure thumb
- Advanced Ozone™ material, promotes air flow
- Breathable in general use, providing good skin care and around-the-clock comfort in different heat and humidity conditions
- Anatomically designed, providing tailored individual fit
- Lightweight, low profile design
- Available in two sizes

Ordering Information

Part Number	Wrist Circumference		Size
	CM	Inches	
82-96191	14 - 18 cm	5½" - 7"	T1
82-96192	18 - 22 cm	7" - 8½"	T2

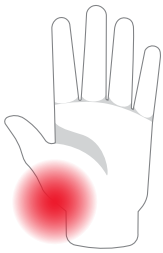

DonJoy CarpaForm®

Functional carpal tunnel wrist brace

Indications


- Carpal tunnel syndrome
- Tendinopathy, Teno-synovitis

Characteristics

- Ergonomic, ultra-lightweight design
- Supportive exoskeleton design for decompression of carpal nerve
- Malleable aluminum structure provides tailored support
- Supports the wrist in neutral position, can be used at night
- Dedicated for either right or left hand use
- Easy and simple to adjust and fit
- Available in three sizes

Ordering Information

Part Number		Wrist Circumference		Size
Right	Left	CM	Inches	
82-96232	82-96222	13 - 16.5 cm	5" - 6½"	S
82-96233	82-96223	16.5 - 19.5 cm	6½" - 8"	M
82-96234	82-96224	19.5 - 25 cm	8" - 10"	L

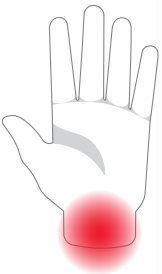

DonJoy RhizoForm®


Functional malleable thumb brace

Indications

- Indicated for rhizarthrosis of the thumb

Characteristics

- Ergonomic universal design, fits both right and left hands
- Ultra-lightweight supportive exoskeleton structure
- Malleable structure provides tailored support
- onSteam® breathable microfiber lining promotes comfort
- Easy and simple to both adjust and fit
- Available in two sizes

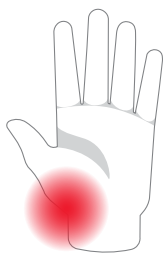


Ordering Information

Part Number		Wrist Circumference		Size
Right	Left	CM	Inches	
82-96292	82-96242	13.5-16.5 cm	5½" - 6½"	S/M
82-96293	82-96243	16.5-20 cm	6½" - 7½"	L/XL



The Technology of Comfort

We realize the importance of the form, fit, and long-term comfort requirements that users have for a brace, both as a medical device and as a wearable product within everyday lifestyle needs. We listened to our users and ensured that this comprehensive range of wrist and hand braces benefited from the use of the most appropriate materials. These are incorporated throughout the range and include:

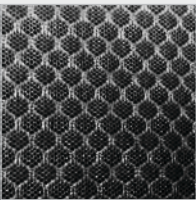


Ozone™

- Honeycomb patterned 3D spacer mesh material, with soft lining
- Innovative and highly breathable, promotes airflow and enhances patient comfort



Flexy™

- Highly breathable 3D spacer material combined with a velvet exterior, for comfort

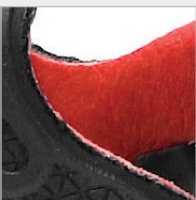

onSteam®

- Microfiber moisture evacuating, thermo-regulating lining
- Antibacterial, anti-odour, and washable, to help maintain hygiene
- Hypoallergenic, skin-friendly properties
- Moldable in conjunction with exoskeleton for improved fit

Note: Not all materials used on all products



DONJOY®

Meet the family.

This range comprises nine innovative braces covering the gamut of hand and wrist injuries, from sprains and strains to tendinopathies and rhizarthrosis. By incorporating valuable customer feedback into our development process, we've been able to maximise the benefits of this highly functional and stylish offering. Available in a variety of sizes and configurations, patients can be confident that no matter their need, we've got it covered.



DUOFORM



RESPIFORM



CARPAFORM

WRIST

WRIST (Carpal tunnel)



DUOFORM+



RESPIFORM+



ERGOFORM



RHIZOFORM

WRIST & THUMB

THUMB



DIGIFORM

FINGERS (Maximum 2)



DIGIFORM+

FINGERS (Maximum 3 or 4)



DONJOY®

DJO Global | 1a Guildford Business Park | Guildford | Surrey | GU2 8XG | UNITED KINGDOM

www.DJOglobal.eu

DONJOY DUOFORM, DONJOY DUOFORM+, DONJOY RESPIFORM, DONJOY RESPIFORM+, DONJOY DIGIFORM, DONJOY DIGIFORM+, DONJOY ERGOFORM, DONJOY CARPAFORM and DONJOY RHIZOFORM are trademarks of DJO, LLC. All other trademarks depicted herein are the property of their respective owners and there is no sponsorship, association, or affiliation between DJO, LLC and those trademark owners.

Individual results may vary. Neither DJO Global, Inc. nor any of its subsidiaries dispense medical advice. The contents of this brochure do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.